

DST-1A

Economy type
Electric integrated operating table



Unique advantages >

The footboard is detachable and easy to adjust, which is very suitable for urology surgery; the tabletop can be moved horizontally and transparently. The tabletop can move horizontally, tilt forward and backward, tilt left and right, and rise and fall are all realized by electric push rods, which are operated by buttons and are safe and reliable.




ADVANTAGES AND FEATURES

- ① Adopting precise medical micro-motor and electric push rod mechanical structure, the product has superior rigidity, uniform and stable running speed, accurate body position, reliable performance and durability.
- ② Lifting, translation, head and foot tilting, left and right tilting, back lifting and other actions are driven by independent power systems.
- ③ The built-in waist bridge is greater than or equal to 120mm.
- ④ The bed surface has a large electric horizontal movement stroke of 300mm, and it can be adjusted appropriately during surgery without moving the patient. High-precision linear guide rails ensure the smooth operation of the operating table.
- ⑤ Mattress: Made of a whole sponge with high sealing and soft texture. Natural shaping, with memory recovery function, reducing local pressure. Seamless, waterproof, easy to clean, anti-static, and removable.
- ⑥ The operating bed frame is made of high-quality nickel-chromium stainless alloy steel. The table side rails are made of stainless steel 304 and then surface sandblasted and electrolytically oxidized.
- ⑦ The base shell is galvanized steel plate sprayed, which is beautiful and easy to clean and disinfect.


SPECIFICATION AND SIZE




Headboard up

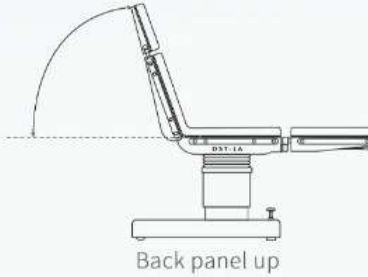


Back panel up

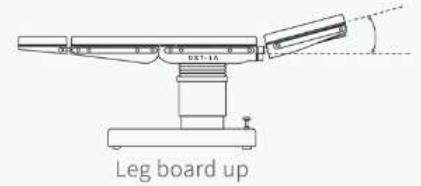


Leg board up

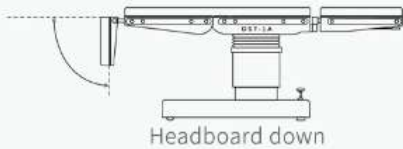


Headboard down

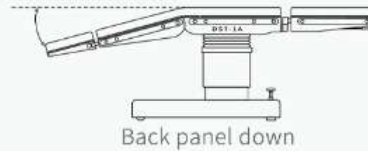


Back panel down

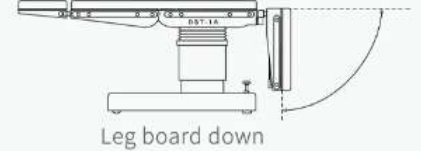


Leg board down

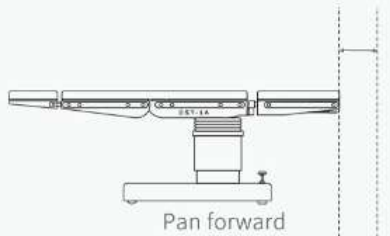


Pan forward

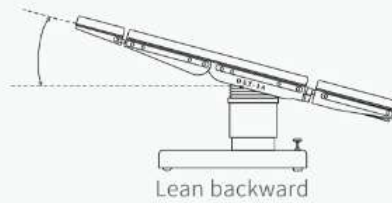


Lean backward

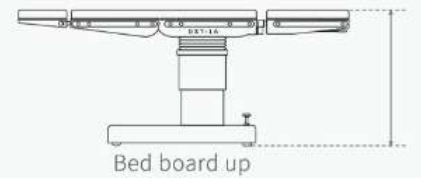


Bed board up

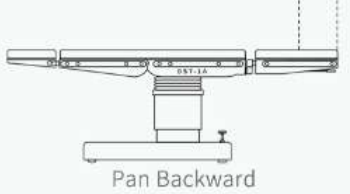


Pan Backward

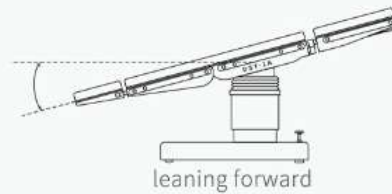



leaning forward

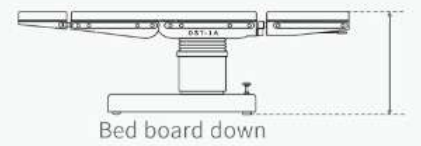

Bed board down

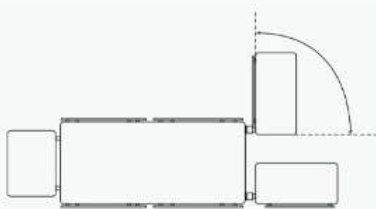

Leg board moves outward

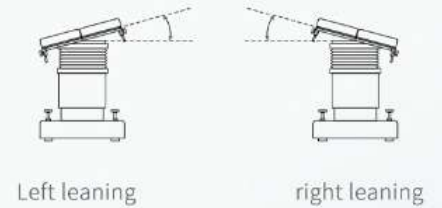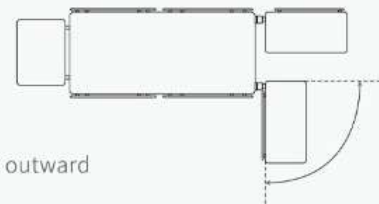

Left leaning

right leaning

STANDARD ACCESSORIES

Anaesthesia screen	1 piece	Body support	1 pair	Leg support	1 pair	Hand remote	1 piece
Kidney bridge handle	1 piece	Arm support	1 pair	Foot support	1 pair	Mattress	1 set
Shoulder support	1 pair	Long fixing clamp	1 pair	Fixing clamp	8 pieces	Power Line	1piece